




# Surshield™ Versatus™-WP

## IV catheter with injection port

### Instructions for use


#### 1 Preparation

- Prepare skin according to local policy and make sure all items are ready.
- Remove the IV catheter from the packaging and remove the cover from the needle just before use.

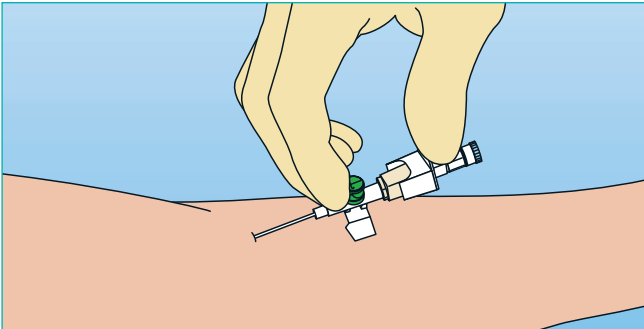


#### 2 Venepuncture

- Take the catheter with the intended transparent housing. **Do not** hold or touch the ivory protection cover at any moment.
- Insert the needle into the vein using low angle and check flashback.
- Lower angle and advance the catheter slightly.

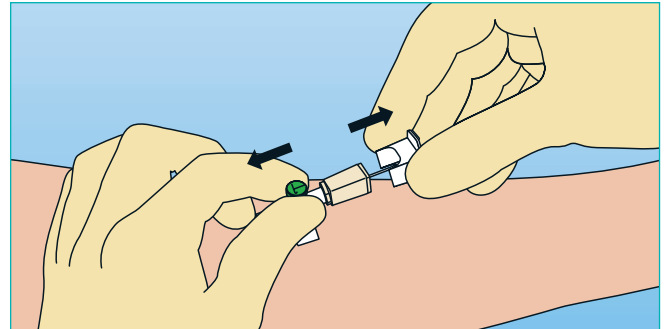



#### 3 Insertion

- Thread the catheter into the vein while removing the needle. Note that either a one handed or two handed technique can be used.

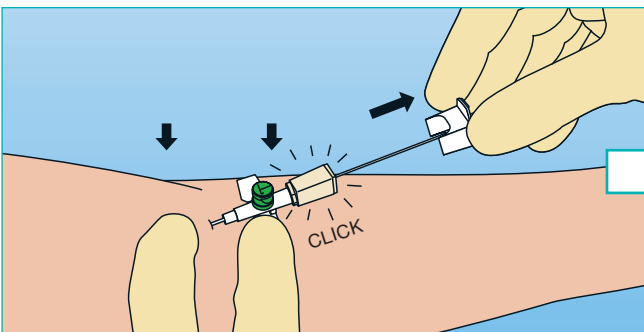

#### 4 Needle removal

- Apply pressure above the catheter tip, stabilize catheter hub and remove the needle from the catheter hub. The needle point is encapsulated automatically into the protection cover. The user is now protected against accidental needlestick injuries.
- **Do not** hold or touch the ivory protection cover at any moment.

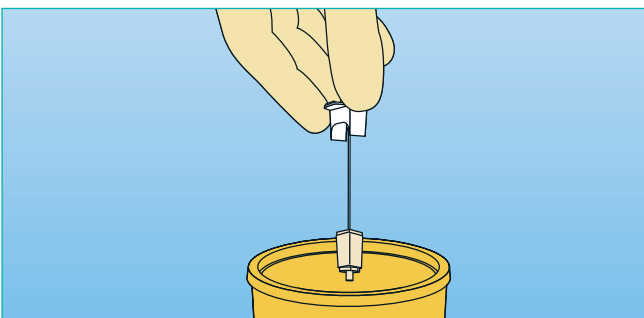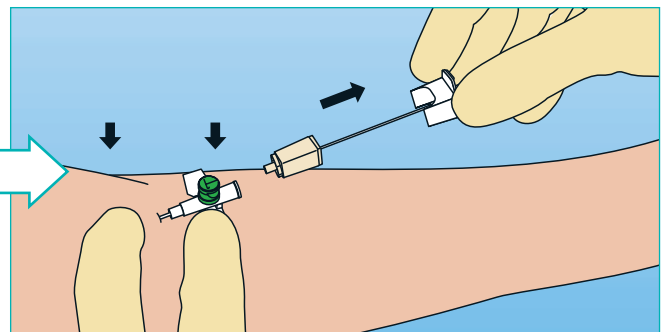

#### 5 Disposal

- Dispose of the protected needle into a sharps container immediately

#### 6 Securement

- Tape and dress per local policy.



We keep life flowing

TERUMO EUROPE N.V. • Researchpark Zone 2 Haasrode  
• Interleuvenlaan 40, B-3001 Leuven BELGIUM  
• Tel. +32/16.38.12.11 • Fax +32/16.40.02.49  
@: Registered Trademark • Published by Terumo Europe N.V.  
www.terumo.com • www.terumo-europe.com




# Surshield™ Versatus™-S/Surshield™ Versatus™-W

## Straight IV catheter/IV catheter with wings

### Instructions for use


#### 1 Preparation

- Prepare skin according to local policy and make sure all items are ready.
- Remove the IV catheter from the packaging and remove the cover from the needle just before use.

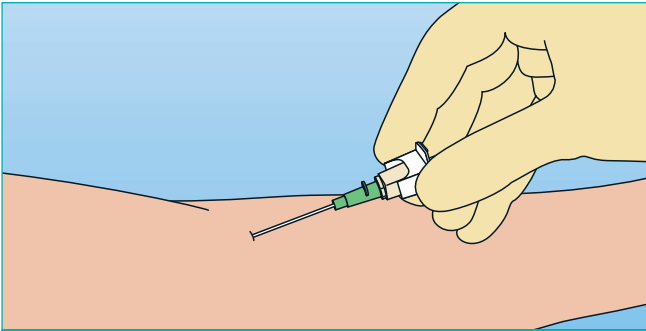


#### 2 Venepuncture

- Take the catheter with the intended transparent housing. **Do not** hold or touch the ivory protection cover at any moment.
- Insert the needle into the vein using low angle and check flashback.
- Lower angle and advance the catheter slightly.

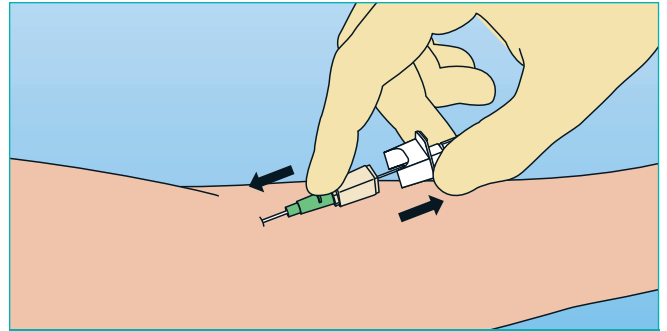



#### 3 Insertion

- Thread the catheter into the vein while removing the needle. Note that either a one handed or two handed technique can be used.

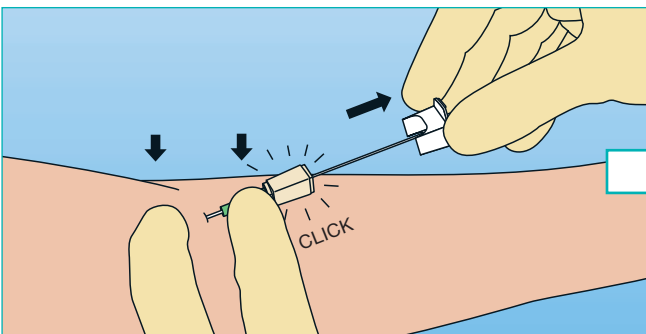

#### 4 Needle removal

- Apply pressure above the catheter tip, stabilize catheter hub and remove the needle from the catheter hub. The needle point is encapsulated automatically into the protection cover. The user is now protected against accidental needlestick injuries.
- **Do not** hold or touch the ivory protection cover at any moment.

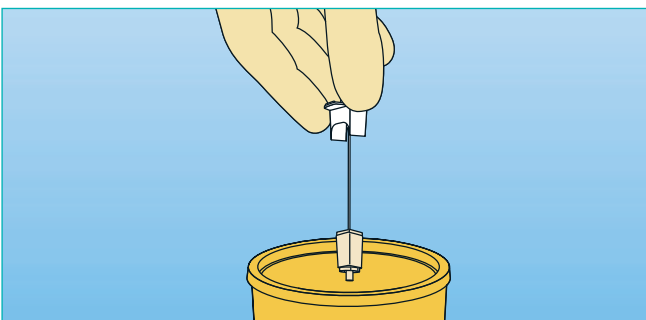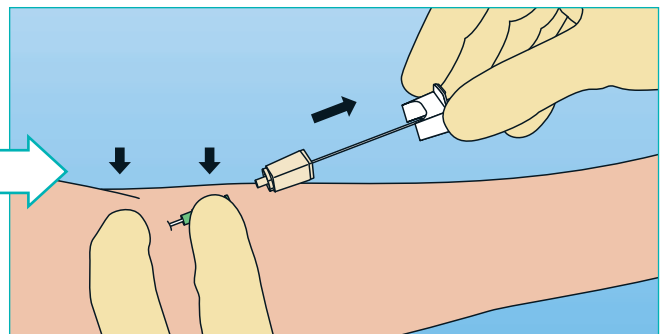

#### 5 Disposal

- Dispose of the protected needle into a sharps container immediately

#### 6 Securement

- Tape and dress per local policy.



We keep life flowing

TERUMO EUROPE N.V. • Researchpark Zone 2 Haasrode  
• Interleuvenlaan 40, B-3001 Leuven BELGIUM  
• Tel. +32/16.38.12.11 • Fax +32/16.40.02.49  
@: Registered Trademark • Published by Terumo Europe N.V.  
www.terumo.com • www.terumo-europe.com